

Test Menu

Thyroid

FT3
FT4
T3
T4
TSH
Anti-TG
Anti-TPO
Tg
rT3
TRAB*

Bone

Metabolic
Intact PTH
Calcitonin
Vitamin D total

Adrenal Function

DHEA-S
Cortisol
ACTH

Hypertension

Angiotension I*
Angiotension II*
Renin*
Aldosterone*

DS' Screening

PAPP-A*
Free β-HCG*

Infectious Disease

HIV Combo
HBsAg
Anti-HBs
HBeAg
Anti-HBe
Anti-HBc
Anti-HCV
Anti-TP (Syphilis)
HBeAg quantitative*
Anti-HBc IgM*

Cardiac

CK-MB
Myoglobin
Troponin I
BNP

Liver Fibrosis

LN
HA
PIIINP
CIV

Growth Hormone

hGH*
IGF-1*

*: in development

Fertility

Total β-HCG
FSH
LH
Prolactin
E2
E3
TESTO
PROG
Free testosterone*
17-OH PROG*
SHBG*
AMH*

Diabetes

Insulin
C-peptide

ToRCH

Toxo IgG*
Toxo IgM*
Rubella IgG*
Rubella IgM*
CMV IgG*
CMV IgM*
HSV-1/2 IgG*
HSV-1/2 IgM*
HSV-1 IgG*
HSV-2 IgG*

Tumor marker

CEA
AFP
CA125
CA15-3
CA19-9
FPSA
TPSA
NSE
CYFRA 21-1
CA72-4
PG I
PG II
SCCA
HE4
ProGRP
CA50
CA242

Anemia

Ferritin
Vitamin B12
Folate
RBC Folate

Inflammation

PCT

Technical Specifications

General Information

Throughput Up to 180 tests/hour
Measurement Principle Micron superparamagnetic particles platform
With alkaline phosphatase (ALP) labeled reagents
And AMPPD substrate

Sample Assembly

Sample Feeder Continuous sample loading and offloading
Sample positions 50
Sample Probe Steel probe, liquid level detection, clot detection,
Horizontal and vertical collision protection
Sample Volume 10~200ul, with increment of 1ul
Sample Dilution Automatic dilution with 1:2~1:40

Reagent Assembly

Reagent Disk 15 reagent positions, continuous loading,
real-time mixing for magnetic particles
Cooling System 2~8 °C
Reagent Pack 50 tests/pack, 100 tests/pack
Reagent Scanning Built-in/External bar code reader
Reagent Probe Steel probe, liquid level detection, clot detection,
Horizontal and vertical collision protection
Reagent Volume 10~200ul, with increment of 1ul

Substrate System

Substrate Bottle 300/500 tests per bottle, bottom piercing
Substrate Loading Capable of continuous loading during testing
Substrate Heating Constant heating
Aspiration Volume 200ul

Magnetic separation unit

Separation Assembly 3-phase magnetic separation
Temperature Accuracy: 37±0.3 °C

Measurement and Reaction System

Detection Mode Photon counting
Signal Detector Photomultiplier (PMT) photo counter
PMT Calibration LED reference module
Incubation positions 82
Temperature Accuracy: 37±0.3 °C
Mixing Unit Non-touch vortexer mixing

Cuvette Loading and Waste Collection

Cuvette Disposable plastic cuvette
Cuvette Loading 88 cuvettes per tray,
2 trays can be loaded in one batch,
continuous loading and offloading
Waste collection Waste container, 200 cuvettes per box

Working conditions

Power Supply 100-240V~50/60Hz
Input Power 500VA
Operating Temperature 15~30 °C
Relative Humidity 35%RH~85%RH, without condensation
Dimension 860mm*740mm*560mm (W*D*H)
Weight 132Kg
Altitude Height -400m~3000m (70kpa)

CL-960i

Chemiluminescence Immunoassay system



mindray





LAB Market

بنها - عمارة ٢ غرب الاستاد - بجوار كنتاكي - خلف نقابة الاطباء ومطعم الرشيدى

لاب ماركت

- مستلزمات معملية
- كيمواويات
- أجهزة
- صيانة

01096411333



LAB Market

بنها - عمارة ٢ غرب الاستاد - بجوار كنتاكي - خلف نقابة الاطباء ومطعم الرشيدى

لاب ماركت
 • مستلزمات معملية
 • كيمائيات
 • اجهزة
 • صيانة

01096411333

One of the world's smallest,
fully automated CLIA analyzers



CL-960i The right size for your lab



Precise magnetic separation with consistent performance reliability




One-key start
Extremely easy to use



Solid-liquid waste separation
Environmental friendly



Intuitive software interface
Easy access to all functions

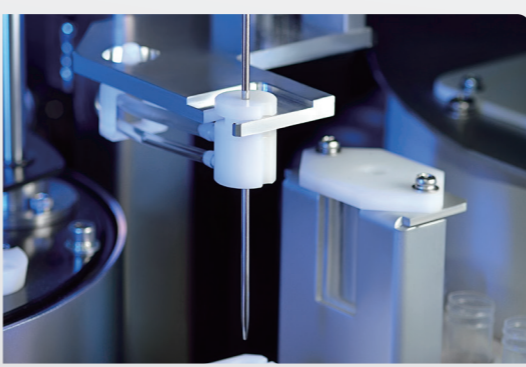


Wiper washing technology
Enhance efficiency and minimize carryover



Continuously loading of reagents and consumables



User maintenance free
Lighten the workload

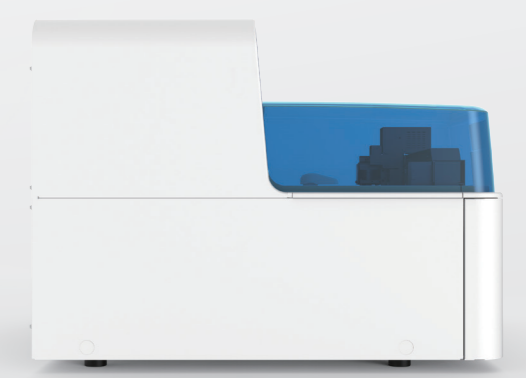

Intelligent consumables management